


PICK IT SMART

SAFE® TUBE PRESENTER

Retrieve your samples faster than ever before and minimize thawing cycles.

- Flexible and durable: Compatible with various rack formats and built to last for daily lab use.
- Simple and safe sample handling: Individual tubes or entire columns can be removed quickly and error-free.
- Efficient workflow: reduces laboratory workload and minimizes the risk of errors.

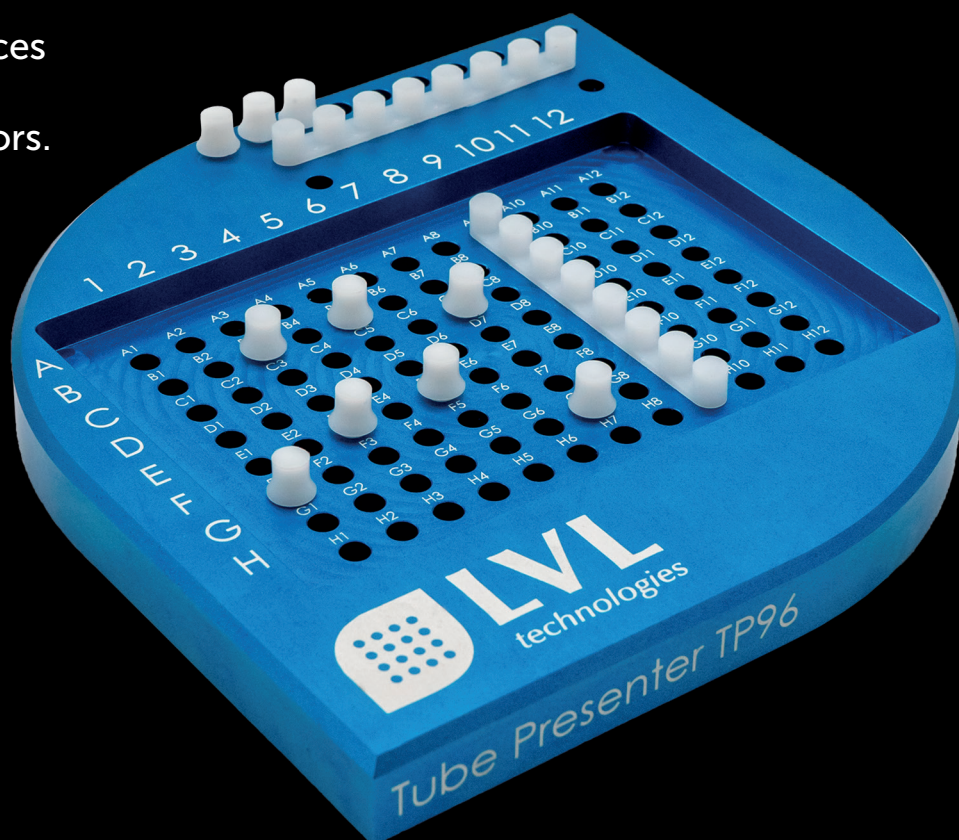

order now.

<https://lvl-technologies.com/en/extras/tube-presenter/>

LVL SAFE® Tube Presenter

The **Tube Presenter** from **SAFE®** is the simple tool for processing and checking out individual samples: simply insert the single pin into the appropriate position, place the rack on top, and the upright tube can be removed safely and easily. The column pin can be used for entire columns, which simplifies the handling of entire sample series. This well-designed solution reduces the workload for laboratory staff and contributes to a smooth workflow by minimizing the risk of errors and maximizing productivity. The **Tube Presenter** is compatible with various rack formats and offers flexible and reliable support for everyday laboratory work. With its robust and durable construction, it is an indispensable tool for any laboratory that values precision and efficiency.



SAFE® Tube Presenter 24



SAFE® Tube Presenter 48



SAFE® Tube Presenter 96



SAFE® Tube Presenter 81



SAFE® Tube Presenter 100

Technical information

Version / Article No.	TP24	TP48	TP96	TP81	TP100
Compatible with following SAFE® racks	SBS24	SBS48	SBS96	Cryo Box 9x9	Cryo Box 10x10
Description	Includes 4 single pins and two 6-pins connectors	Includes 4 single pins and two 6-pins connectors	Includes 10 single pins and two 8-pins connectors	incl. 20 single pins and two 9-pins connectors	incl. 20 single pins and two 10-pins connectors
Dimension L x W x H	160 x 160 x 15 mm	160 x 160 x 15 mm	160 x 160 x 15 mm	220 x 200 x 15 mm	220 x 200 x 15 mm